

Deliver Reliable Results with Maximum Uptime

H60 Hematology Analyzer



DETAILS THAT LIGHTEN YOUR WORK



Reagent management via RFID card secure your reagent business



12-inch colorful touch screen




Two routine reagents to support test



16 microliter of sample aspiration

INTUITIVE SYSTEM THAT SIMPLIFIES YOUR HANDLING

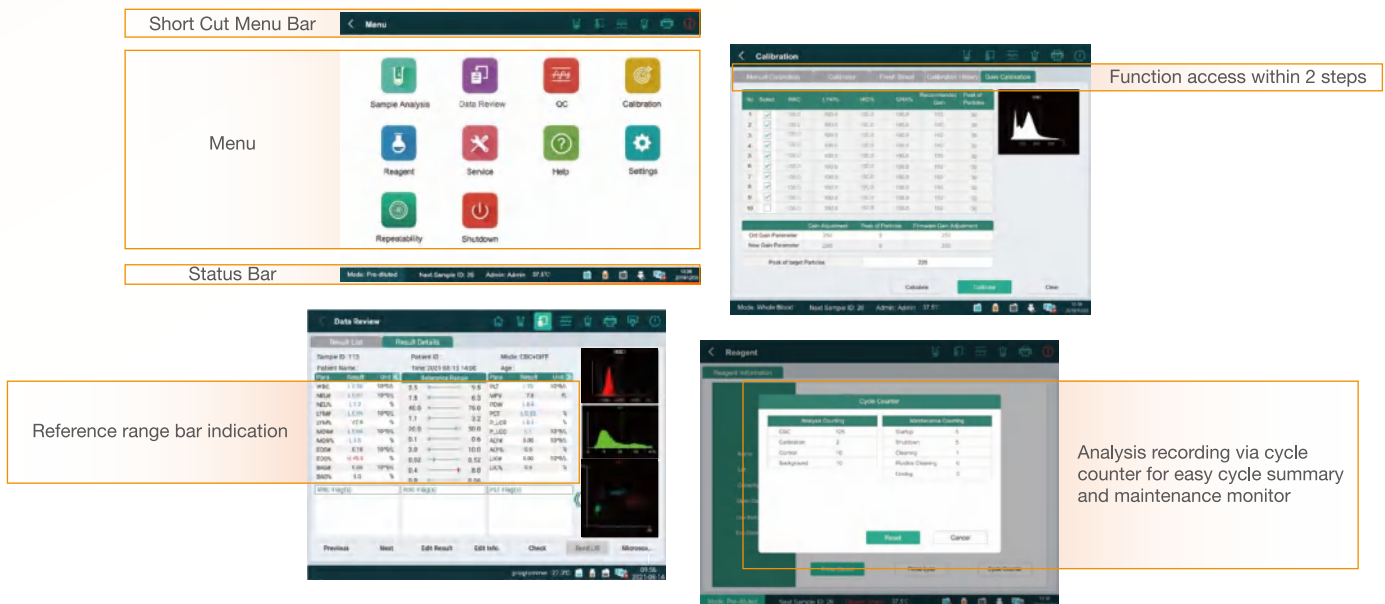

The system that stands out

The H60 features a whole set of user-oriented functions that are hard to find elsewhere. Its android-based system has a similar experience as a smartphone – pleasure-to-operate and easy-to-use. You can get to any functionality within two steps. The 12-inch colorful touch screen makes working with H60S and viewing its results exceptionally straightforward.

Outstanding graphic interface

The user interface is undoubtedly the most visible and best-appreciated feature of the H60 series. Different set-up views, each optimized for various user situations, facilitate routine use. The most common include a sample view showing classic graphs and all parameters within one page. The other sample view features the most practical reference bar display that indicates the results' status. The comprehensive information on the short-cut menu bar and status bar give you fast access to any primary function.

Communicate and print like never before


An integrated workstation with the HIS and LIS connections further streamlines analysis and reporting. Standardized HL7 protocol and storage capacity extends up to 100,000 easily accessible sample files. The print function has the most flexible options that allow you to print individual sample results and control tests, L-J graphs, or QC lists with any parameters selection.

SENSITIVE FLAGGINGS THAT ASSIST YOUR INSIGHTS

NRBC

NRBC PRESENCE FLAGGING

REMINDED TO CORRECT WHITE BLOOD CELL COUNTS BY SMEAR

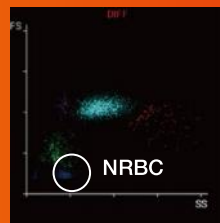

ALY

ABNORMAL LYMPHOCYTE FLAGGING


SUSPICIOUS OF VARIANT LYMPHOCYTES AND MALIGNANT LYMPHOCYTES



EOSINOPHILIA

HIGH SPECIFICITY ON EOSINOPHIL


GREAT DIFFERENTIATION ON EOSINOPHILIA SAMPLES

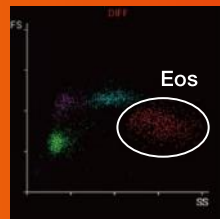

IG

IG PRESENCE FLAGGING
ASSESSMENT OF SEPSIS AND INFLAMMATION

COMPRISED OF METAMYELOCYTES, MYELOCYTES & PROMYELOCYTES

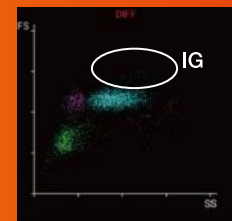


DIFF


10 DIFF PARAMETERS

HIGH CLUSTERING OF 5-DIFF

TWO 2D SCATTERGRAMS
HIGH SENSITIVITY FLAGGING

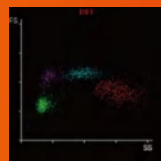

DIFF

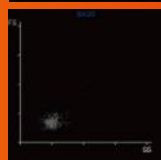

BASO



Specification

Parameters NEW		
25 reportable parameters:	WBC, Lym%, Mon%, Neu%, Bas%, Eos%, Lym#, Mon#, Neu#, Eos#, Bas#, RBC, HGB HCT, MCV, MCH, MCHC, RDW-CV, RDW-SD, PLT, MPV, PDW, PCT, P-LCR, P-LCC.	
8 research parameters:	LIC%, LIC#, ALY%, ALY#, NRBC% , NRBC# , NLR,PLR	
3 histograms:	WBC , RBC, PLT	
1 3D scattergram:	DIFF	
4 2D scattergram:	DIFF, BASO	
Performance		
Parameter	Linearity Range	Precision (CV%)
WBC (10 ⁹ /L)	0.0-500.0	≤3.0% (3.50-7.00) ≤2.5% (7.01-15.00)
RBC (10 ¹² /L)	0.0-8.5	≤1.5% (3.5-6.5)
HGB (g/L)	0-250	≤1.5% (100-180)
PLT (10 ⁹ /L)	0-4000	≤5.0% (100-500)
MCV (fL)		≤1.0% (70-120)
Principles		Throughput
Tri-angle laser flow cytometry analysis for WBC, DIFF and BASO counting Electrical impedance method for RBC, PLT counting Cyanide-free reagent for HGB with colorimetric method		60 samples per hour
Sample Volume		Interface
Whole blood mode	16 μL	12.1 inch colorful touch screen
Capillary whole blood mode	16 μL	
Pre-diluted mode	20 μL	
Reagent		Control and Calibrator
HD600 Diluent 20L		ED-60D, ED-CAL PLUS
HL600 Lyse 500mL/1L		Data Storage Capacity
HC310 Cleaner 50mL		100, 000 results including results and histograms 60 QC files (100 data per file)
Operating Environment		Dimension and Weight
Temperature: 15°C~32°C; Humidity: 30% RH~85% RH; Air pressure: 70 kPa~106 kPa		501(D)×320(W)×502.5(H) Weight: 29.3kg

Hospisolution Limited
Tel.: + 86 13912963835

Web: www.hospisolution.com
Email: info@hospisolution.com