

YA-PS06B

Hydraulic Stretcher



*Better care
Better comfort!*

YA-PS06B

Hydraulic Stretcher

DESCRIPTION

YA-PS06B is a hydraulic patient transfer stretcher for inside hospital patient transferring. The hydraulic operation, safety guard rail, corner bumpers, central brake system etc.. all the featured functions make the stretcher an ideal option for hospital and medical institutions at an affordable budget.

- FEATURES
- Height adjustable in a hydraulic way by pedals on both sides
 - Collapsible aluminum alloy railing
 - Central brake system
 - Pushing bar located on both ends
 - Back section adjustable in a pneumatic way
 - Bumper castor located on four corners of the stretcher

STANDARD SPECIFICATION	Overall size	2110*690*460mm
	Back section adjustment	0-75°
	Height adjustment	460-810mm
	Caster	Diameter 125mm
	Locking	Side control central locking
	Loading Capacity	180kg


Structure

All the base and frame structure is made by thickness steel. With dust cover on the bottom, it is durable and easy to clean and disinfect.


Foot pedal

Hydraulic pump controlled by foot pedal with anti-slip cover.


Mesh bed

2-section mesh bed platform with good air permeability.


Castor

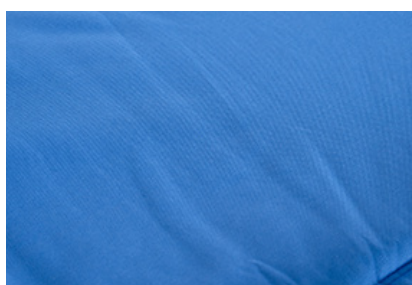
Dia 125mm castor with central locking, side control for easy operation.


Side rail

The 5-bar long side rail prevents the patient from rolling over while moving.


Push handle

Ergonomics push handle in head and foot end position.


Mattress

Mattress with 4 handles for easy strapping or patient transfer.


Bump castor

Oxygen tank holder

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US